




AJL Capsular tension ring



Model	AC / ACS				
References	AC001002	AC001102	AC001202	AC001302	AC001402
	ACS001002	ACS001102	ACS001202	ACS001302	ACS001402
Open Diameter	10mm	11mm	12mm	13mm	14mm
Close Diameter	9mm	10mm	10mm	11mm	12mm
Brand name	AJL CAPSULAR TENSION RING				
Product description	AJL capsular tension rings are PMMA medical devices used during cataract surgery to stabilize the capsule. The rings have different diameters and sections (AC/ ACS) to accommodate patients' individual capsular bags.				
Indications	Indicated in cases of weak zonulas where there is risk of the capsule is retracted. The purpose of these implants (AC/ ACS) is to expand the capsule in which is located the crystalline to allow the introduction of an intraocular lens inside. Suitable for implantation together with an intraocular lens in aphakia after the extraction of the crystalline. The square section of ACS ring prevents the migration of cells into the posterior chamber reducing posterior chamber opacification (PCO).				
Material	PMMA (Polimethylmethacrylate) Medical Grade with UV filter				
Design	AC: round section / ACS: square section				
Presentation	Each ring is placed on a holder specially designed to fit each model to enable the ring extraction for its implementation. The holder avoids any damage during transportation. Not preloaded.				
Sterilization Method	Ethylene Oxide (ETO)				
Expiration	4,5 years				

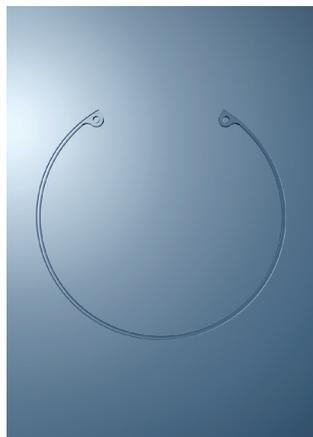